

Jazzy[®]

SELECT

#1 IN EASE OF USE.
#1 IN MANEUVERABILITY.

The highly maneuverable Jazzy[®] Select[™] delivers the advantages of in-line motor technology for enhanced efficiency, torque, range, and performance. Plus, it's loaded with performance and convenience features, like standard Active-Trac[®] Suspension, that make it as easy to enjoy as it is to use.

#1 Physician Approved

JAZZY[®] SELECT[™]
shown in Jazzy Red with
gray vinyl contemporary design,
high-back seat, PG GC controller
and foot platform



Authentic
Jazzy

Pride[®]
Mobility Products Corp.

800-800-8586 (US) • 888-570-1113 (Canada) • www.pridemobility.com

JZ SELECT 09/07



The Jazzy® Select™ series models include:

- Jazzy Select-C
- Jazzy Select-SS

FEATURES

- Dual in-line motors for enhanced efficiency, torque, range, and performance
- New 40 amp, PG GC controller
- Built with ease of service in mind
- Easy side access to batteries
- Tool-less front panel access to electronics to minimize service time


OPTIONS

- High-back seat with headrest
- High-back seat with solid seat pan
- Cup holder
- Weather cover
- Rear basket
- Oxygen tank holder
- Cane/crutch holder
- Walker holder

The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed and range vary with user weight, terrain type, battery charge, battery condition and tire pressure. Due to manufacturing tolerances and continual product improvement, length, width, turning radius and ground clearance can be subject to variance of + or - 3 %.

Jazzy® Select™	
Weight capacity:	300 lbs.
Maximum speed:	Up to 4 mph
Ground clearance:	2"
Turning radius:	20"
Overall length:	35.5" without front riggings
Overall width:	23.75"
Wheels:	Drive: 10"
Front anti-tips:	5" solid
Rear casters:	6" solid
Suspension:	Active-Trac®
Drive train:	Two-motor, in-line, mid-wheel drive
Braking system:	Intelligent braking (electronic, regenerative disc brakes)
Std. Electronics:	40 amp, PG GC controller
Specialty controls:	N/A
Battery charger:	O -board, 3 amp
Per-charge range:	Up to 15 miles
Battery requirements:	(2) 12 volt, deep cycle
Battery size:	U-1
Transport weight	
Battery weight:	24.5 lbs. each
Base weight:	99 lbs.
Max. Seat dimensions:	20" x 20"
Warranty:	5-year limited warranty on frame; 1-year limited warranty on electronics; 1-year limited warranty on drive motors

Colors:

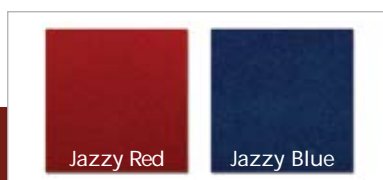

Pride SilverStar Lifts & Ramps

- Backpacker® Plus
- Full Platform
- Power chair lift with straps



Backpacker® Plus



Full Platform

These lifts are compatible with the Jazzy® Select™.